

MK MAYAKOBA

**Model: KAN-TCHR-2307**  
**SHIATSULOGIC DX Massage Chair**

**Operational manual**

*Make your life healthy  
and beautiful every day!*



Please read this operating instruction carefully before using this chair  
Please keep this manual in perfect condition

## Please keep this manual in perfect condition

Thank you for purchasing our massage chair. To ensure optimal function and safety, please read this manual carefully.

## Content

Important safety instruction .....	2
Preventative safety measures .....	4
Parts name and function .....	5
Function Introduction .....	5
Operation Method .....	6
Parts name and function of remote control .....	7
Operation prevention measures .....	9
Grounding instruction .....	9
Cleaning and maintenance .....	10
Trouble shooting Guide .....	11
Specifications .....	12

## Important safety instruction

For of this item, make sure to follow basic prevention measures including the following points:

### Environment for operating the chair

- Please make sure to leave 40cm of space the chair, so the backrest can recline.
- Do not operate this item in direct sunlight or in a place with high temperature, exposure could fade or harden leather.
- Please do not operate this item in damp areas such as bathroom to prevent electrical shock or other accidents.
- Please keep this item safely away from heat sources, humidity, sharp objects, corrosive chemicals, and flammable objects.
- Please do not use chair outdoor.

### Before using this item

#### Turn on the power after checking the surrounding.

1. Check the surrounding:

- Check if there are any objects behind the chair.
- Check if there is enough space (at least 40cm) so the chair can be reclined.

2. Check item:

- Check if there is a cleft on head pad, central back pad. also check cloth material.
- Check if the supply cord is broken or if there is dust on the plug.

3. Turn on the power.

- Make sure the voltage conforms with the required specification of this item.
- Insert the plug into the power socket.  
turn on the power switch on the side of the chair backrest.

## Precaution measure during use of the chair

**Dangerous!** To reduce electric shock risk to minimum:

1. Remember to unplug the power after each use.
2. Remember to unplug the power before cleaning.

**Warning:** to reduce the risk of burning damage, fire, electricity shock, or bodily harm:

1. Before using the chair, make sure to uncover the back pad to check if there is any tear or holes in the cloth. If any tear or holes found, stop using immediately and send it to authorized repair center.
2. The plug of this item can only be connected to a suitable grounded socket. Please refer to grounding instructions. The plug should be accessible after installation.
3. This item should only be operated according to this manual.
4. Please don't stand on the chair.
5. Please remember to unplug the chair if not in use. Please close all the massage functions before unplugging the chair.
6. If the supplied cord is damaged, it must be replaced by the manufacturer or an authorized service agent to avoid damage or injury.
7. In the chair becomes damaged, splashed or dropped into water, or fails to operate, please contact the dealer and send it to an authorized service center for repair.
8. To avoid electric shock or bodily harm, please don't open or disassemble the chair by yourself.

## The people who are not suitable for using the chair

1. Disabled people or those who can't express themselves must be strictly and continuously supervised while using this chair.
2. Please don't allow children to use or play around the chair.

**Warning:** to avoid bodily harm:  
Preventative safety measures

1. Pregnant women, people of poor health and those who feel pain in their neck, back, waist, or hip should consult doctors before use. People with heart pacemakers or other personal medical devices should consult doctors before use.
2. Don't put any part of your body between the rollers, otherwise roller's massage action may cause harm. Also do not put your hands or feet into the gap between the back and the seat, back and armrest or other active parts to avoid bodily injury.
3. Do not use massage mechanism against bare skin. Although you may receive a stronger massage effect, it may cause skin irritation. Do not wear any hard objects while massage.
4. Do not massage the head, belly, elbow or knees. Please be cautious when massaging shoulder areas.
5. Do check the position of your neck, if the position of neck is not suitable, please adjust to the suitable position for safety.
6. Maximum massage time is 15 minutes (automatic stop function after 15 minutes). A massage of more than 15 minutes may cause an uncomfortable feeling.
7. If you feel nausea or pain during the massage, please stop the massage, you may press "stop" key at any time, it will stop massage immediately.
8. Please don't sleep while massaging. Massage after drinking is forbidden.
9. When back of the chair is in a reclined position, please don't sit on it or press heavily on it, as it may flip upside down.
10. Don't allow children or pets to play on it especially during massage.
11. If power problems occur, please unplug the power immediately.

## Parts name and function

Foot spa massage chair



## Function Introduction

1. Equipped with ultra thin 2D massage machine on back, which can achieve kneading, knocking, shiatsu, tapping and manipulation etc. massage methods to simulate a professional massage therapist.
2. With the body shape detection functions, it is suitable for people in different height and size to carry out an accurate massage.
3. Combined with ergonomics curved guide rail, you can fully enjoy the head to thighs massage, to quickly relieve the spinal pressure.
4. Vibration massage at seat can help relieve fatigue of thighs.
5. Armrests can be lifted up/down for easy access, flexible handle on right armrest for hanging bag, very convenient.
6. Humanization design LCD controller offers a real time display of massage position, methods and status.
7. With mobile phone charging function.

## Operation Method & Instruction

Below are some guidance for using this chair. For further details, please refer to the operation instruction of remote control.

1. Switch on.
2. Press on/off key and select massage technique.
3. With the body shape detection functions, it is suitable for people in different height and size to carry out an accurate massage.
4. Press other keys on the remote control to choose massage functions that you need.  
When you think the massage strength is too weak, you may remove the pillow or back cushion to get a stronger massage.

Automatic Massage:

- Very suitable for those whole body feels tired.
- Very suitable for those who need simple automatic massage.

Manual massage

- To choose some specific massage function or fixed massage area.

5. Back and Seat Adjustment:

- Back Adjustment: press the button of back-up or back-down to the ideal position, then release it; The back-rest's moving angle is 110°~135°.
- Seat Adjustment: press the button of move-forward or move-backward to the ideal position, then release it; The seat moving distance is 0~130mm.

6. Vibration massage at seat can help relieve fatigue of thighs.

8. With mobile phone charging function.

9. Switch off and massage is over.

- Press ON/OFF key to turn off all massage functions and wait for massage roller restored.
- Switch off and unplug the power cord.


During the process of whole massage, if the pressure on the massage roller is too large, for security, this machine will start the internal safety protection circuit immediately to stop all the massage action. The corresponding indicator light of the controller will flash continuously, the buzzer will produce "di-di" constant warning sound, the controller button will reject all the commands.

- Indicator light of "Kneading" is flashing: means the kneading mechanism is overloaded.
- Indicator light of "Tapping" is flashing: means the tapping mechanism is overloaded.
- Indicator light of "Overall" is flashing: means the up & down mechanism is overloaded.

If the above problems happen, please turn off the power switch and wait for 10 seconds, restart the power switch and do the operation. If you feel the massage pressure is too hard or it moves slowly, please raise your body a little (this can happen frequently for people over 100kg, please use the machine carefully to avoid permanent damage and hurt on your body).

# Parts name and function of remote control

Operation instruction of remote control:

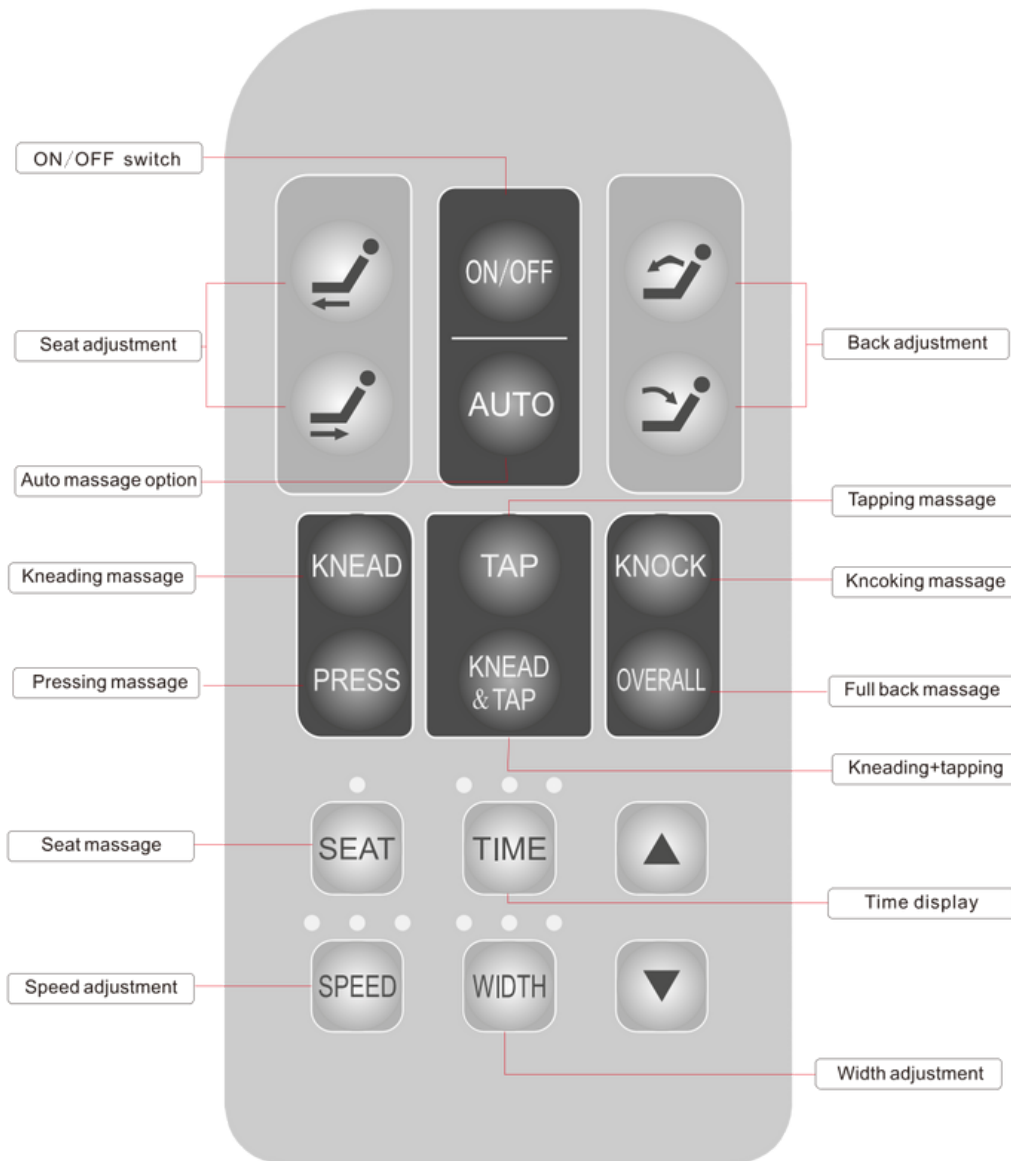

- ON/OFF** Press this button to switch on /off in turn. When power on, Automatic icon and Massage Method icon on the LCD screen will be lit up; when power off, STOP icon on the LCD screen will be lit up.
- AUTO** Press this button, the massage machine enters or exits the automatic control state. There are three modes under the automatic state, each time press the button, followed by the back massage - shoulder massage - waist massage - automatically stop, the corresponding image on LCD screen is flashing. In auto mode, seat pressure or vibration massage function is opened automatically.
- TIME** The default display is 20mins, each time press the time key, the time display is followed by 10 → 20 → 30mins. After time is set ok, the display time decrements every minute until to 00, then the machine will shut down.
- KNEAD** In manual state, you can manually open and close a variety of massage techniques. Press this button, massage machine starts kneading, the corresponding light-emitting diode lights on and the kneading image is flashing on the LCD screen.
- TAP** Press this button, massage machine starts tapping, the corresponding light-emitting diode lights on and the tapping image is flashing on the LCD screen.
- KNEAD & TAP** Press this button, the massage turns into knead + tapp mode, the corresponding diode lights up and the image of knead + tapp on the LED screen is flashing.
- PRESS** Press this button, the massage turns into pressing mode, the corresponding diode lights up and the image of pressing on the LED screen is flashing.
- KNOCK** Press this button, the massage turns into knocking mode, the corresponding diode lights up and the image of knocking on the LED screen is flashing.
- OVERALL** Press this button, the massage machine enters the whole work state, the corresponding light-emitting diode is bright; The corresponding image on the LCD screen is flashing. (In automatic state, this key is invalid.)
- SPEED** In manual mode, the working state of the massage machine is kneading, tapping, knead&Tap, knocking and shiatsu, press this key to adjust the speed by turns, 1→2→3→4→5. The corresponding display speed bar on the LCD screen is gradually increasing.
- WIDTH** In manual mode, if the working state of the massage machine is tapping, shiatsu and knocking, press this button, width can be adjusted to narrow→middle→wide by turns. The corresponding image is displayed on the LCD screen.
- Seat adjustment icons** Keep pressing this button for a while, seat slides forward; release the key, seat stops moving.
- Seat adjustment icons** Keep pressing this button for a while, seat slides backward; release the key, seat stops moving.
- Back adjustment icons** Keep pressing this button for a while, back-rest raises up; release the key, back-rest stops moving.
- Back adjustment icons** Keep pressing this button for a while, back-rest moves down; release the key, back-rest stops moving.
- SEAT** Press this button, the corresponding indicator light will be on, vibration function or air massage on the seat will start to work and the corresponding seat icons on the LCD screen starts to flash.
- Up arrow** In manual state, press this button for a long time, massage machine will go up; Press this button for a long time, massage machine will go down.
- Down arrow** In automatic state, these two buttons are invalid.

## Operation prevention measures

Warning: to avoid damage to the chair

1. User's maximum weight is 120kgs; weight excess may cause damage to massage machine or permanent damage to other parts.
2. Please don't crush or drop the control set. Don't pull control cable or electric wire. Be careful not to clamp the cable between armrest and back of the chair.
3. Please don't touch cloth material of the chair with sharp things, and be careful not to drop cigarette, ash or lighters on the chair. Chair exposure in the direct sunlight may fade the cloth material.
4. If there is abnormal function in the chair, please shut off power and ask the designated repair center to check the chair. Don't repair it by yourself.
5. Don't repeat pressing any keys in very short time in case of any abnormal functions.
6. Noise coming from inside when it is turned on is normal; it's the item's internal structure.
7. Make sure there are no obstacles behind the chair before adjusting the chair's position. The maximum adjust range is 0-130mm. Please lower the chair slowly and gradually to avoid touching the wall or other things that can make the chair function improperly.
8. Each time when finish the massage, please turn off the power and unplug the chair.

## Grounding instruction

- To reduce the risk of electric shock, proper grounding is required for the chair.
- This item is equipped with a grounded plug and special grounding bolt; insert the grounding plug into a grounded socket with proper installation and suitable local voltage, or connect the special bolt to reliable objects with grounding wire.

### Dangerous!

- Incorrect grounding may cause electrical shock. If you have any questions about this item, please consult a licensed electrician. If the plug sent along with the chair doesn't match your socket, please don't modify the socket. Ask technicians to install a suitable one. Modification from the original manufactured specifications will void all warranties expressed or implied.
- Please don't use a plug adapter.

## Cleaning and maintenance

### Synthetic leather

- Clean the leather with a piece of soft dry cloth.
- If the leather becomes dirty, please mix water with neutral 3%-5% detergent and dip the cloth into the water for cleaning the leather surface. Then wash the cloth to remove detergent and scrub the leather with a piece of dry cloth. Let it dry naturally (don't dry the wet leather with an electric dryer.)
- If leather has been contacting ethylene resin for a long time it might fade.

### Cloth cover

1. Use a piece of cloth with medium detergent liquid to clean cloth cover.
2. Use a soft brush dipped with medium detergent to brush dirty cloth material. Be sure not to use too much, it might cause damage to the cloth.
3. Dip a piece of soft cloth into clean water. Use it to wash the item and let it dry naturally.

### Cleaning of plastics parts

1. Dip the cloth into medium detergent liquid.
  2. Use a piece of cloth with clean water to scrub this item and dry it naturally.
- For cleaning of the remote control, please use a slightly damp cloth to scrub the surface.
  - Don't use any solvents to clean as they may cause discoloration.

## Caution:

- Be sure to unplug the chair when cleaning this item.
- After cleaning this item, make sure to dry hands before plugging the power again. Otherwise, it may cause electric shock or burn.

## Trouble shooting Guide

<p>This item can't be operated at all.</p>	<ul style="list-style-type: none"> <li>● Electric wire has been taken off.</li> <li>● Power key is still on OFF mode.</li> <li>● Massage roller is on the way of storage .It will work after the roller comes to the storage position.</li> <li>● Electric wire or plug is broken. It should be replaced by manufacturer, appointed repairman, or qualified technicians.</li> <li>● Fuse is broken. Please replace it with a new one of same specification.</li> <li>● For internal circuit troubles, please send it to appointed repair place for check.</li> </ul>
<p>Massage roller stops working while in operation;control indicator keeps flashing;buzzer gives out continuousalarm.</p> <ul style="list-style-type: none"> <li>● <b>Kneading indicator flashes:</b> Relative operating mechanism of kneading has problem.</li> <li>● <b>Tapping indicator flashes:</b> Relative operating mechanism of tapping has problem.</li> <li>● <b>The entire course indicator flashes:</b> Relative operating mechanism of up/down has problem.</li> </ul>	<ul style="list-style-type: none"> <li>● Over pressure and overload on the massage roller make internal circuit start to protect itself, please reduce load and shut off the power to restart the item in 10 seconds.</li> <li>● Internal massage mechanism is damaged. circuit may start its self-protection. Please askappointed department to check.</li> <li>● Broken cloth wraps or prevents massage roller from moving, so the protection circuit starts its protection function. Please ask appointed department to check.</li> </ul>
<p>Too much noise</p>	<p>Long term use may cause mechanism abrasion. it's normal.Please ask the appinted department to check.</p>
<p>Control mechanism becomes too hot</p>	<p>Extended use of control may generate excessive heat.If overheating occurs please send the unit to an authorized service center</p>
<p>Electric wire or plug is abnormally hot</p>	<p>Immediately stop use and send it to the authorized repair facility.</p>

## Specifications

- Model:KAN-TCHR-2307
- Power supply: 120V~ 60Hz
- Power consumption:100W
- Massage area on back:
  - Up/down:about 42cm
  - Left/right:about 19cm
- Backrest decline angle:approx. 110~135°



AYC Group LLC  
 Email: [service@aycllc.com](mailto:service@aycllc.com)  
[www.aycllc.com](http://www.aycllc.com)